

Case Report


Olecranon Fracture

Melvin P. Rosenwasser M.D.

Columbia University New York

August 2006

- Pre-Op

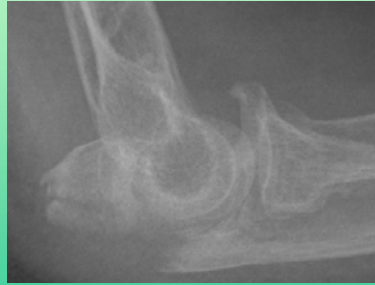


Displaced olecranon fracture
(Mayo Type IIA)

Case Report

Olecranon Fracture

- Step 1:


Insert 2mm K-wires
into ulna and reduce
fracture

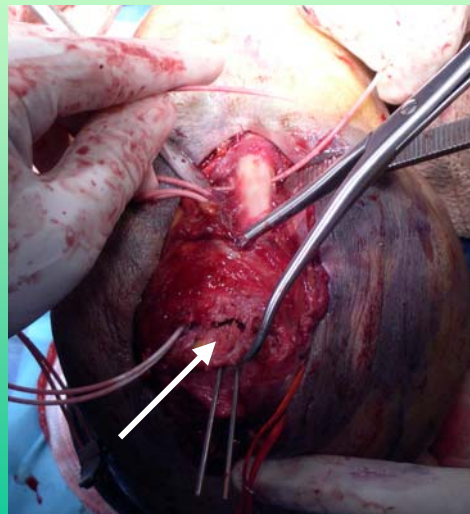

Case Report
Olecranon Fracture

- Step 2:


Drill transverse
hole through ulna
and pass cable

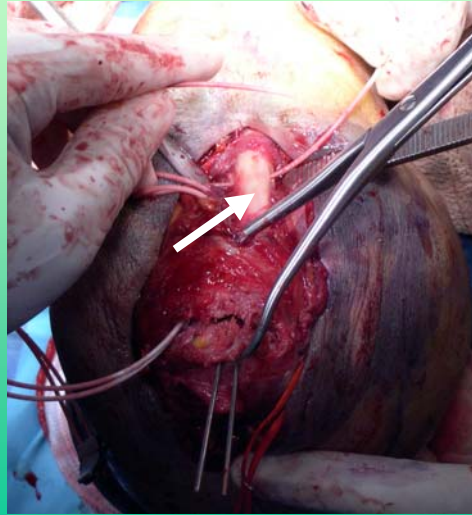

Case Report
Olecranon Fracture

- Step 3:

Wrap cable around K-wires

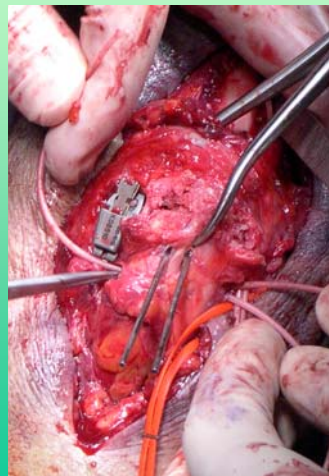


Case Report
Olecranon Fracture

- Step 4:

Position clasp under
triceps on lateral side


Apply tension and lock
cables

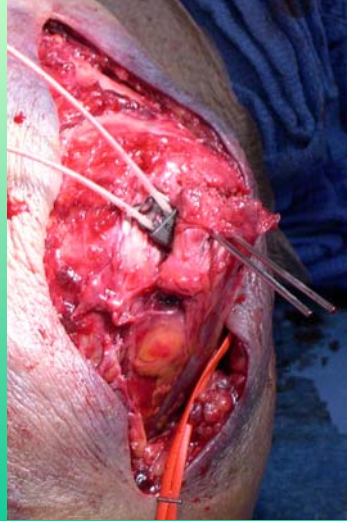

Case Report
Olecranon Fracture

- Step 5:

Triceps closure

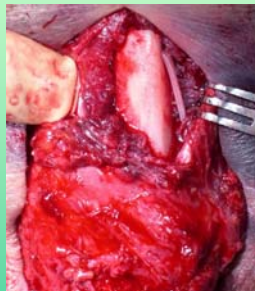

Case Report
Olecranon Fracture

- Post-Op

Hardware retained

Patient healed
uneventfully

